

AEL-100

MEDICAL EXAM LIGHT IN ALLTION QUALITY


ALLTION AEL-100 LED MEDICAL EXAM LIGHT
HIGHER FIDELITY, MORE FLEXIBILITY AND POSSIBILITIES STATIONARY AND
YET FLEXIBLE, BRIGHT AND YET HOMOGENEOUS






BRIGHTER, COST SAVINGS AND LESS WASTE

- 1 white 3W LED, typical light output and pinpoint accuracy.
- Life time up to 50,000 hours. High Illuminance 70,000 Lux at 200MM working distance.
- Less energy consumption.
- * No bulb replacement.

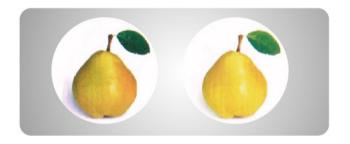

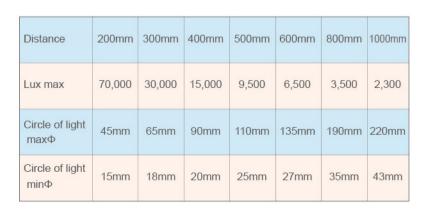



Conventional LED ALLTION LED

AUTHENTIC COLOR RENDERING

- The light is more pleasant to the eye.
- Excellent color rendering. The illminated area appears natural and contrasted.



ADJUSTABLE SPOT SIZE

- Spot size can be adjusted in a diameter between 15mm to 220mm to suit any situation at wide-range working distance 200mm to 1000mm.
- · Easy to clean.



FLEXIBLE GOOSE NECK

- Extremely flexible goose neck can be fixed in the chosen position and stopped accurately without spring back.
- *Two-segment goose neck design could be bent radius to any angle and in all directions.









ONE-HAND OPERATION HANDLE PANEL

- Akimbo-style handle highlights at a glance, convenient to operate whatever pushing, pulling or rotating only by one hand.
- The prominent brightness adjustment knob and touchable black On/Off are humanized designed and assembled on the handle panel.

A VARIETY OF MOUNTING OPTIONSAVAILABLE

- Mobile stand. Easily to be moved due to 5 low friction lockable universal castors.
- Wall mount.Quickly to be mounted on conventional walls.
- * Table Mount.Simply to be mounted through universal clamp.
- Pole Mount.Flexible clamping at any suitable space circular pole, adjustable clamping diameter (Φ 25-Φ40 mm)













Pole Mount

DC POWER SUPPLY AND EASY ASSEMBLY

- Optional DC power supply.With DC5Vpower adapter (Providing US,European,British, and other international plugs).Can use reliable DC5V power bank for supplying.
- The prominent brightness adjustment knob and touchable black On/Off are humanized designed and assembled on the handle panel.

Technical Parameters

Light Specifications	LED	1 White 3W LED
	Life Time	50,000 hours
	Color Temperature	5,300K
	Spot Diameter Adjustable @ Working Distance 200mm	15-45mm
	Illuminance @ Working Distance 200mm	70,000lux
Physical Diameters	Goose Neck Length	1000mm
	Stand Pole Height	700mm
	Base Diameter	500mm
	Gross Weight	6KGS
	Net Weight	3.5KGS
	Package Measurement	86×61×16(cm)
Electrical	Voltage	DC 5V
	Power	5W
	Power Cable	5.5×2.1mm
	Adapter	Input : AC100-240V ~ 50Hz/Hz
		Output : DC 5V
Misc Data	Mounting Options	Mobile Stand, Table / Wall / Pole Mount
	Kind of Extension	Goose Neck
	Warranty	2 Years
	Usage Environment	5°C-40°C , 30%-80%RH,860hpa-1060hpa
	Storage Environment	-5°C-40°C , 30%-80%RH,860hpa-1060hpa

Alltion reserve the right to change specification without notice.



MICROSCOPES ARE OUR BUSINESS & SERVICE IS OUR PRIDE

Alltion (Guangxi) Instrument Co.,Ltd.

Alltion (Wuzhou) Co.,Ltd.

Alltion Building, NO.10,3RD.Road, Wuzhou Industrial Park, Wuzhou, Guangxi, China 543000

Tel: 0774-2836101,2836102

Fax: 0774-2836192

Email:sales@alltion.com or sales@alltion-microscope.com Web site:www.alltion.com or www.alltion-microscope.com