ISO13485



AC-2000 SERIES COLPOSCOPE

FEATURES AND BENEFITS

Apochromatic Optic Design
3-Step Magnification and/or Single Magnification
LED Illumination for Long Life & Low Cost
1080P FULL HD Real-time output module
Two Optical Circles for Quick Measuring Size of Tissue
Integrated Green Filter
Compact, Light Weight Base with Wheels for Easily Maneuver
Ergonomic Design for Long Time Precision Operation Comfortably
Completely Assembled and Ready to Use






SUPERB OPTICAL FEATURES

Our excellent optical design makes the colposcope with superb optical feature features, including:

- 1. High resolution
- 2.Real 3-Dimensional Precise Image Reproduction
- 3.Large depth of field
- 4. Wide Field of View
- 5.Good Contrast
- 6.1080P FULL HD Real-time output of image

Apochromatic optics design could correct the chromatic aberration of three types of colored light effectively. It ensures the reduction and saturation of image color and improves the resolution, contrast and view depth of image.

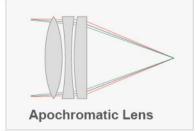

Color Aberrations from Visible Lights Are Corrected



Visible Light Refraction froms Color Aberrations


3-STEP MAGNIFICATION AND/ OR SINGLE MAGNIFICATION

3-step magnification (3.75x-7.5x-15x) provides larger magnifying range and more selection on magnification to meet the requirement of your practice, from lowest magnification for an overview to high magnification for precise and close inspection and treatment.

For matching your budget requirement, we supply single magnification model as well.



TWO OPTICAL CIRCLES FOR QUICK MEASURING SIZE OF TISSUE

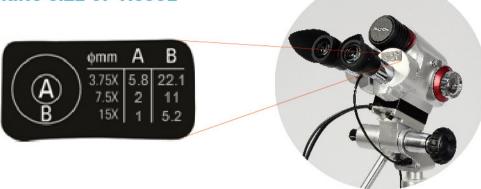

INTEGRATED GREEN FILTER

Green filter can help you to recognize blood vessels.





LED ILLUMINATION

LED coaxial illuminator provides pure white color light similar to the sun light to show the real color. The bulb life can be over 10,000 hours/10 years for lasting performance. It not only avoids the trouble of changing bulbs, but also saves your cost for the spare bulbs & electric charge.

INTEGRATED FULL HD VIDEO CAMERA

An integrated 1080P FULL HD video camera is available in specific models. It is very convenient for the doctors to communicate with the patient by this device. The images from integrated video camera can be inputted into a monitor or a TV set directly. Or it can be inputted into a computer by image collection with USB adapter or stored into a hard disk recorder.





COMPACT, LIGHT WEIGHT BASE WITH WHEEL FOR EASILY MANEUVER

The colposcope is light to be lifted easily only 11kgs. The base comes with 5 wheels so that it stands stably and moves freely.

For matching your working space, we supply tilt base model as well.





ERGONOMIC DESIGN FOR LONG TIME PRECISION OPERATION COMFORTABLY

All ergonomic design assure you can work with our colposcope for long time precision operation most comfortably. It would avoid or reduce your problems on the back and neck.

COMPLETELY ASSEMBLED AND READY TO USE

Alltion AC-2000 series colposcopes are delivered to customers with completely assembled unit in one carton. Any installation is not necessary. All you have to do is plug it in and/or connect it to your monitor or /TV set.



SPECIFICATION ALLTION AC-2000 SERIES COLPOSCOPES:

- 1.Total magnification 3.75x, 7.5x, 15x for 3-step models and 7.5x for single magnification models
- 2.300mm focal length
- 3.Field of View: 79mm, 39mm, 19mm dia. for 3-step magnification models and 39mm for single magnification models
- 4.Depth of Field: 4.5mm, 1.13mm, 0.76mm for 3-step magnification models and 1.13mm for single magnification models
- 5.Light intensity: >30,000 LUX
- 6.Individually adjustable eyepiece 16.7x
- 7.Right eyepiece with built-in two circles for you to measure the size of the problem you see.
- 8.Built-in rheostat for brightness adjustment
- 9.1080P FULL HD Real-time output module
- 10.User-selectable built-in green filter for enhanced contrast
- 11. LED bulb is easily changeable.
- 12. Working Height: 950mm to 1250mm
- 13. Fine focus adjustment Handles
- 14. Two (Gross/Fine) height adjustment Handles
- 15. Unit is completely assembled and ready for use once power is on
- 16. Main Power Supply: 100V to 240V

Item and specification subject to change, please contact Alltion for current production information.

	MODEL NO.	SPECIFICATION		
		Magnification	Binocular Head	Built in Camera
Complete Unit	AC-2110	Single	Straight	Nil
	AC-2111	Single	Straight	Full HD Video Camera
	AC-2310	3-Step	Straight	Nil
	AC-2311	3-Step	Straight	Full HD Video Camera
Assembly Kits	AC-2110 Assembly Kit including colposcope head, head holder, power supply without post and base			
	AC-2111 Assembly Kit including colposcope head, head holder, power supply without post and base			
	AC-2310 Assembly Kit including colposcope head, head holder, power supply without post and base			
	AC-2311 Assembly Kit including colposcope head, head holder, power supply without post and base			
Complete Unit	AC-2110 Tilting star	Single	Straight	Nil
	AC-2310 Tilting star	3-Step	Straight	Nil
	AC-2311 Tilting star	nd 3-Step	Straight	Full HD Video Camera





MICROSCOPES ARE OUR BUSINESS & SERVICE IS OUR PRIDE

Alltion (Guangxi) Instrument Co.,Ltd

Alltion (Wuzhou) Co.,Ltd

Alltion Building, No. 10, 3rd Road, Wuzhou Industrial Park,

Wuzhou, Guangxi, China. 543100

Tel: +86-0774-2836101,2836102

Fax: +86-0774-2836192

Email: sales@alltion.com or sales@alltion-microscope.com Web site: www.alltion.com or www.alltion-microscope.com